

protect Epico ROM[®]

Post-Op Range of Motion Elbow Brace

Indications

- conservative care of elbow instabilities, with or without other injuries
- orthotic care for medial or lateral epicondyle fractures
- postoperative care after prosthesis implantation
- postoperative care after elbow stabilization
- functional follow-up treatment for stable fractures
- passive mobilization after elbow arthrolysis
- postoperative care of elbow ligament instabilities

Mode of action

- immobilization and stabilization of elbow joint
- flexion and extension limitation can be set in steps of 10° from 0° to 120°
- movement restriction orthosis for early functional treatment

Features

- open and flat frame construction
- total weight of only 400g
- adjustable length handle
- malleable frame
- comfortable shoulder strap included
- latex Free breathable open cell foam padding
- easy to adjust QuickSet hinge
- safety cover protects hinge ROM settings
- low profile flat hinge
- 30% lighter than previous version



protect Epico ROM

Features



30 % less weight
The orthosis has an open and flat frame construction and a total weight of only 400 g.



Adjustable hand cushion
Telescopic hand support to ensure the orthosis is securely attached to the arm.



Comfortable shoulder strap
The added shoulder strap ensures the wearing comfort for the patient.




Latex free padding
The orthosis has latex free breathable open cell foam padding.



Low profile flat Quick Set hinge
The quick set hinge is easy to adjust and has safety covers to protect the hinge ROM settings.

protect Epico ROM

size	length cm	article number left	article number right
standard	up to 29	99450	99350
long	29 -42 cm	99451	99351

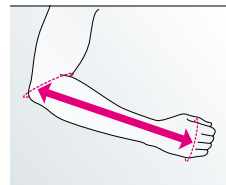

measurement
measure length from elbow to knuckles